



HIGHLIGHTS

- The Live Virtual Collaboration Kit allows a KARL STORZ NEO IP integrated operating room to add bi-directional audio, video and Keyboard/Video/Mouse Switching (KVM) control of an owner-furnished computer running Teams/Zoom/Webex or any other cloud-based virtual collaboration platform. Surgical teams can easily communicate from the OR with anyone in the world virtually and focus on what matters during the patient's procedure.
- Plug-and-play addition that uses standard IT conferencing applications already deployed on your secured hospital network, requiring very little effort from IT departments to deploy.
- Removes specialized/redundant conferencing applications from the operative department, instead using standardized and cost-effective IT solutions.
- Any video source in the operating room may be sent to the far-end viewer in single, PIP, side-by-side or quadview modes.*
- The PC, account licensing and any additional security features of the cloud conferencing platform are owner-furnished.

CONNECTIVITY

Owner-Furnished Equipment (O.F.E.) Computer Requirements	1x HDMI or DisplayPort Graphics Output 1x 3.5 mm Stereo Audio Output Quad Core 3.2 GHz or better 8 GB RAM (16 recommended, depending on how many other apps the computer will run) 2x USB 3.0 or better with no other high-bandwidth USB devices on the USB 3.0 controller chipset Windows 7, 10, 11
Features	1x Audio/Video Feed into communications platform 1x Audio/Video Feed from computer into OR integration system Telestration, chat, emojis, etc., as provided by owner's virtual conferencing platform of choice.


COMPATIBILITY

Integration Systems	OR1 KARL STORZ NEO® IP v2.2.2, v2.3.1, v2.4+
Cloud Communication Platforms	Agnostic to cloud communication platform

* PIP, side-by-side, and quadview may require the purchase of additional integration features depending on your current OR configuration.

It is recommended to check the suitability of the product for the intended procedure prior to use.
Please note that the described products in this medium may not be available yet in all countries due to different regulatory requirements.

DIAGRAM

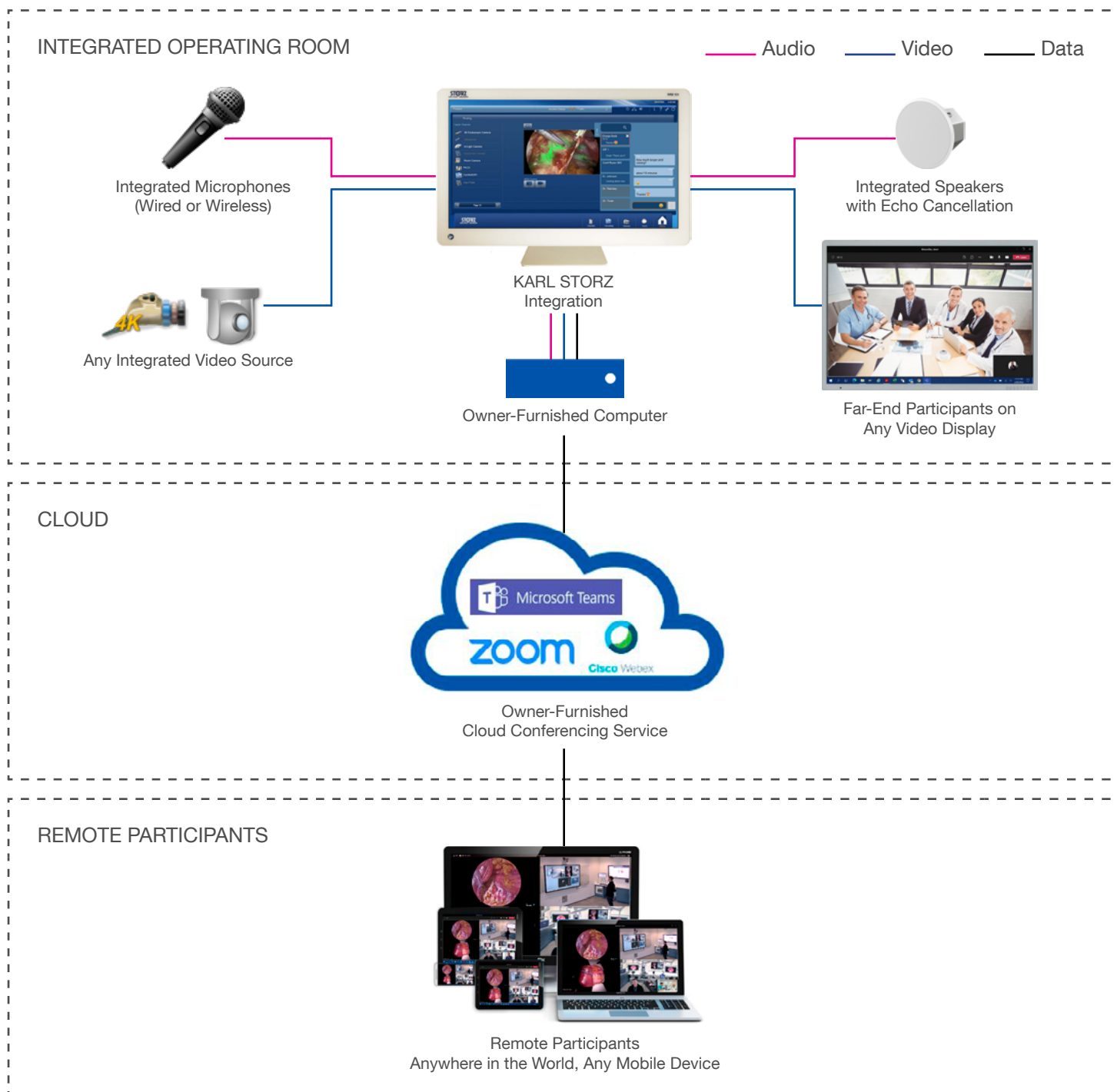

It is recommended to check the suitability of the product for the intended procedure prior to use.
Please note that the described products in this medium may not be available yet in all countries due to different regulatory requirements.



KARL STORZ Endoscopy-America, Inc.
2151 East Grand Avenue | El Segundo, CA 90245-5017, USA
Phone: 424-218-8100 Toll-free: 800-421-0837
Fax: 424-218-8525 Toll-free: 800-321-1304
E-Mail: communications@karlstorz.com

www.karlstorz.com

FOLLOW KARL STORZ

